



ZTM:A1™ Transport and Processing Media (Supplemented with Antibiotics: Clindamycin and Gentamycin)



Product Names: ZTM:A1™

Product Codes: ZTMA1-100, ZTMA1-500, ZTMA1-500-6

General Use: Media to transport and process tissues and cells

Features: GMP quality; sterile; USP grade materials; no animal-derived components; no human proteins; with antibiotics

Manufacturer: INCELL Corporation LLC
Address: 12734 Cimarron Path
City/State/Zip: San Antonio, TX 78249
Phone: 800.364.1765; 210.877.0100
Email: info@incell.com
FAX: 210.877.0200

Distributor: Charter Medical
Address: 3948-A Westpoint Boulevard
City/State/Zip: Winston-Salem, NC 27103
Phone: 866.458.3116 (Toll Free)
Email: incell@chartermedical.com
FAX: 336.714.4241

General Description

ZTM:A1™ transport and processing media is a sterile, animal component free, chemically defined, GMP-manufactured, and ready-to-use solution to ship and process freshly isolated tissues and separated cells or tissues intended for clinical use. ZTM:A1™ is a proprietary formulations of USP and pharmaceutical grade chemicals, with selected antibiotics, and supplemented with a long-acting glutamine derivative to support maintenance and viability of tissues during transport or processing. ZTM:A1™ media does not contain supplemental growth-promoting factors, or human or animal-derived sera or supplements, so it is not suitable for cell culture.

Formulation and Packaging

ZTM:A1™ media contains components (salts, amino acids, sugars) that are high quality USP or pharmaceutical grade chemicals from plant and microbial sources. ZTM:A1™ media helps maintain living tissue during transport, processing and short term storage. ZTM:A1™ antibiotics are: clindamycin (6.5 µg/mL); and gentamicin (50 µg/mL). There is no amphotericin B, in contrast to ZTM:A2™ which has twice ("2X") clindamycin (13 µg/mL); gentamicin (100 µg/mL) and also contains amphotericin B (2.5 µg/mL). The ZTM:A1™ product is packaged in 100 mL (ZTMA1-100) and 500 mL (ZTMA1-500) sizes in square PETG media bottles that meet medical device quality for media storage.

Use and Methods: Transportation and Processing

ZTM™ media help maintain a supportive tissue and "cell-friendly" environment during transport and subsequent bioprocessing (e.g., enzymatic and mechanical tissue disruption). Antibiotics are included to control bio-burden from source tissues or that may have been introduced into the sample at harvest. Table 1 shows tissues where ZTM™ media are suitable for use as the ZTM:A1™ or ZTM:A2™ formulas. For processing tissues, and depending on the protocol, ZTM:A1™ and ZTM:A2™ can be used from 2°C to 37°C. ZTM™ solutions are NOT cryopreservation or culture media. Tissues that are to be processed for living cells, or the released cells, should NOT be frozen or cultured in these media.

Manufacturing

Media are manufactured by sterile 0.22 µm filtration and packaging, using cGMP standards in an ISO Class 7 clean room and ISO Class 5 biosafety cabinet. Raw materials are pre-tested and the final product is checked by quality specifications and acceptance criteria (see below) tested by USP standards for pH, osmolality, sterility (bacteria, fungi), mycoplasma and endotoxin prior to release and preparation of the Certificate of Analysis.

Specifications

Visual
pH (USP <791>)
Osmolality (USP <785>)
Sterility: SC (USP <71>)
Sterility: fTG (USP <71>)
Mycoplasma (USP <63>)
Endotoxin (USP <85>)
Expiration Date

Acceptance Criteria

Colorless, clear, solution
6.5 to 7.0
290 to 350 mmol/kg
No microbial growth
No microbial growth
None detected
<0.1 EU/mL
18 months

Storage

ZTM:A1™ can be stored at room temperature (15°C to 30°C) for up to a week, but store longer-term refrigerated at 2°C to 8°C. Shelf-life is 18 months from date of manufacture.

ZTM:A1™
Transport and Processing Solution

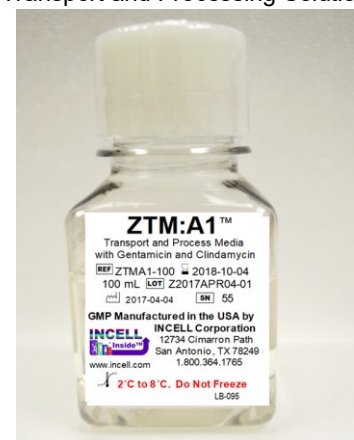


Table 1. ZTM:A1™ and its Related Media* ZTM:A2™ are Suitable for Transport and/or Processing of Surgically Resected Human and Animal Tissues for Use in Various Clinical or Research Applications

Tissue Types	Primary Cell Type(s) Derived
Adipose (Fat)	Mesenchymal Stem Cells; Stromal vascular fraction regenerative cells; adipose cells
Bone Marrow*	Hematopoietic and mesenchymal stem cells; various types of renewable progenitor cells; Endothelial cells; entire population; bone and cartilage progenitors; other cells
Bone	
Brain and Neural (Spinal)	Progenitors; Induction of differentiation
Breast	Primary epithelial cells and/or mesenchymal support cells
Cartilage	Chondrocytes and Progenitors
Colon; Gastrointestinal, incl. oral; esophageal	Primary epithelial and/or mesenchymal support cells or complex tissues in organ-like cultures
Dental Pulp	Hematopoietic and mesenchymal stem cells; Various types of renewable progenitor cells; like bone marrow
Heart	Primary muscle and mesenchymal cells
Kidney	Primary epithelial cells and/or mesenchymal support cells
Liver	Primary epithelial cells and/or mesenchymal support cells
Lung	Primary epithelial cells and/or mesenchymal support cells
Lymph Node*	Lymphocytes and mesenchymal support cells
Muscles (Peripheral; Heart; Smooth)	Pericytes; Mesenchymal or Stromal Stem Cells; other regenerative cells
Pancreas and Other Neuroendocrine Organs	Pancreatic islet beta and acinar cells; other organs (e.g., adrenal)
Parathyroid*	Adenoma cells isolated and stored for potential autologous re-transplantation
Peripheral or apheresis Blood*	Various types of white blood cells, including lymphocytes, macrophages, etc.; Circulating or mesenchymal cells; Endothelial cells.
Placenta	Trophoblasts; Syncytiotrophoblasts; Endothelial cells; Hematopoietic and mesenchymal stem cells; Various renewable progenitor cells
Prostate	Primary epithelial cells and/or mesenchymal support cells
Skin (adult; newborn)	Epidermal keratinocytes; biopsies; hair follicles Dermal Fibroblasts; Mesenchymal cells
Spleen*	Lymphocytes and mesenchymal support cells
Tumors	Epithelial, mesenchymal, lymphoid
Umbilical cord*	Hematopoietic and mesenchymal stem cells; various types of renewable progenitor cells; Endothelial cells

Choice of ZTM™ media type is based on risk assessment of the donor tissue source. ZTM:A2™ has twice the concentration of antibiotics than ZTM:A1™ and also has amphotericin B (anti-fungal). Choice of ZTM:A2™ depends on amount of bioburden and if fungal organisms might be an expected risk. Also, antifungals should not be used with blood derived cells.

*Asterisks next to tissue type recommend ZTM:A1™; remaining tissue types can be either ZTM:A1™ or ZTM:A2™.

Master Files Applications Note

ZTM™ media (ZTM:A2™ and ZTM:A1™) are in FDA Drug and Device Master Files but the media has not been tested by INCELL for any specific diagnostic or therapeutic use. To request use of a Master File call, fax, or send an email to masterfiles@incell.com.

Animal Component Free. INCELL certifies that these media have no “animal-derived components” per the following criteria:

- No animal-derived components are added by INCELL or come from raw materials supplied as components of the Product.
- Product does not come into contact with animal-derived material during manufacturing, processing, handling, or packaging.
- Our Products are manufactured on dedicated animal free equipment in animal free clean rooms.
- This certification applies only to the condition of the above-described Product in its unopened package, and INCELL assumes no responsibility for a Product failing to meet this Statement after handling or use after opening the package.
- Signed “Animal Origin Statement” can be provided on request.

Ordering: Contact Charter Medical

Toll Free: 866.458.3116
incell@chartermedical.com
FAX: 336.714.4241

Technical Assistance

The scientists at INCELL are available to discuss the media or special needs of your cells, and to assist you in your cell culture applications. Call 1-800-364-1765 or e-mail info@incell.com.